

Certification



ISO13485 / MDSAP



Health Canada



CE

Specification

Output Channel	<ul style="list-style-type: none"> Ⓐ 12 Channels consist by two groups; Ⓑ Each channel can be set and adjusted independently.
Treatment Mode	Nine Modes Standard, Toning, Fast Slim, Lymph Drainage, Relaxation, Abdominals, Thighs&Hips, Facial, Bust
Output Intensity	Adjustable, 0-120mA
Treatment Time	1-99min, stepping 1min
Speed (Training frequency)	Slow, Medium and Fast
Safety Protection	Over-current Protection Over-voltage Protection Short Circuit Protection



LGT-2320BE

Electrotherapy Device for Medical Aesthetics



Guangzhou Longest Medical Technology Co., Ltd.

Address: 301 & 401 of Building 2 & Building 3, No.96, Chuangqiang Road, Ningxi Street, Zengcheng District, Guangzhou, Guangdong Province, 511399, China

Tel: 86-020-66353999

Fax: 86-020-66353920

E-mail: export@longest.cn

Website: www.longestmedical.com



Health Life Longest Care

Longest



- ③ Lose weight
- ③ Tone your body
- ③ Get rid of cellulite


Electronic Muscle Stimulation for Beauty

(Model Number:LGT-2320BE) applies low-frequency currents (modulation or non-modulation) to stimulate the motor nerves or muscles to make them contract. After the muscle contractions, the muscle fiber thickens, the weight of the muscle volume increases, the capillaries in the muscle becomes rich, which strengthens the entire muscle to achieve the effect of losing weight and shaping.



Features

- Multiple Treatment Areas
- A 45-minutes session can burn up to 1000 calories
- No weights during your workout
- Intensity automatic locking for your safety
- Non-Invasive, No Downtime

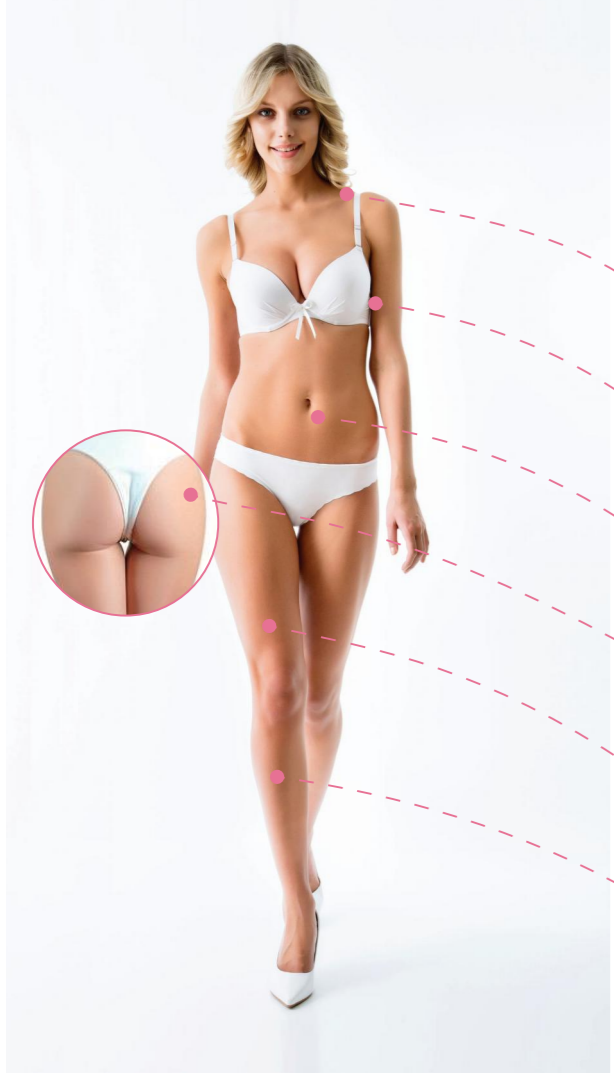

Application

- Back of the Chest

Relaxation Mode
Build up Muscle, Strengthening and Relaxing
- Upper Arms

Standard Mode
Build up Muscle, Strengthening and Toning
- Abdomen

Abdominals Mode
Cellulite Reduction, Strengthening and Toning
- Buttocks

Thigh&Hip Mode
Fat Breakdown and Toning
- Thigh

Thigh&Hip Mode
Break down fat, Strengthening and Lymphatic Drainage
- Lower Leg

Toning Mode
Toning and Cellulite Reduction

Accessory

- 

Silicone Rubber Electrode*12
- 

Straps with 3 Sizes*12
(4 pcs per size)
- 

Electrode Lead Hose
(Φ95mm) *12
- 

Safety Switch*1

